



easeVR

# Getting Started Guide

Brought to you by

appliedVR™



# TABLE OF CONTENTS

This guide will give you a step-by-step overview of EaseVR:

- 1 Daily Practice
- 2 Setting Up the VR Phone
- 3 Setting Up the VR Headset
- 4 How to Play VR
- 5 Using the Heartrate Sensor
- 6 Using the iPhone App
- 7 Starting Your Daily Practice
- 8 Getting Relief
- 9 Troubleshooting
- 10 Turning Off & Charging

## WHAT'S IN YOUR KIT

1



Samsung VR Phone  
and Charger

*Plays the VR sessions*

2



Gear VR Headset

*Allows you to view the  
VR sessions*

3



Heartrate Sensor

*Monitors your heart  
rate to change the  
VR sessions*

4



Headphones  
with Microphone

*Monitors your  
breathing*

## WHAT YOU NEED



Your iPhone

*To use the EaseVR  
app*



EaseVR iPhone App

*Guides you through the  
VR sessions*



Download it by searching for  
"EaseVR" in the App Store

# WELCOME!

*EaseVR uses virtual reality to help you learn essential skills to take control of your pain and provide relief from your pain symptoms. These fun, easy-to-understand and interactive lessons will help you learn how to take manageable steps to change your body and mind's reaction to pain.*

## How do I use EaseVR?

Over the next 8 days, you'll play an EaseVR session each day. Daily practice will help you EASE your symptoms and learn skill you can apply to your everyday life. Daily practice will also provide lasting relief.

## How do I start?

It's easy. Use the EaseVR app on your iPhone to control the VR sessions and follow the 8 chapters of the program. The iPhone app will instruct you on what to play and your VR phone will play the content.

## Why do I need an iPhone and a VR phone?

The app on your iPhone will help you play the VR content and keep you on track. The VR phone is only used for playing VR content.

Learn more at [www.appliedvr.io/ease](http://www.appliedvr.io/ease)

# STEP 2: SET UP THE VIRTUAL REALITY PHONE

## 1. Turn on VR Phone

Turn on the Samsung VR phone by holding down the button on the right side until it vibrates and powers on.



*Power Button*

*Home Button*


## 2. Connect to WIFI

Select the gear icon at the bottom of the screen to open Settings. Tap on Connections, then WI-FI to find your local network. Press the Home button to return to your Home screen.



*Settings*

*Connect to WIFI*

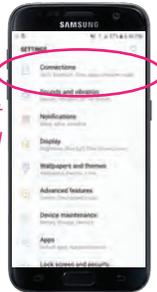

## 3. Add headphones

Plug in the headphones to the bottom of the VR phone.



*Plug in*

# STEP 3: SET UP THE VIRTUAL REALITY HEADSET

## 4. Find the Device Holder

Find the Device Holder on the front of the headset and make sure the latch on the right is pulled back.



*Device Holder with plug*

*Latch without plug*


## 5. Adjust fit

Try on the headset without the phone to make sure the straps fit your head. Tighten side straps first, then tighten top strap.



## 6. Adjust focus

Find the focus dial and volume buttons so you know how to use them when you have the headset on.

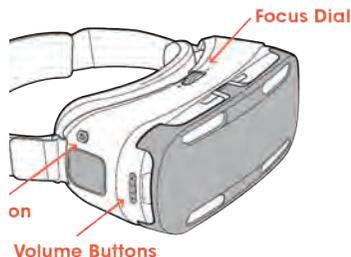

## VR Phone Not Fitting?

If Device Holder is not in B position, push down on the entire highlighted area and slide to B.



*Device Holder in B position*



*Device Holder should stay locked*

## STEP 1: YOUR DAILY SCHEDULE

- Day:**
- 1 Introduction to VR: Observe the moment
  - 2 Start the story & focus on breathing
  - 3 Learn about diaphragmatic breathing
  - 4 Try breathing with a pacer
  - 5 Tree Practice
  - 6 Crystal Practice
  - 7 Black & White Practice
  - 8 Night Fire Practice


**Once you're done with the 8 day program, please feel free to repeat any of the chapters or try the experiences in Relief.**

## STEP 4: PLAY VR

### 7. Open PlayVR

Tap on the PlayVR app on the Home screen.

*PlayVR app*

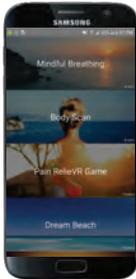

*When you see this screen, you are ready to put the VR phone in the headset*

### 8. Insert VR Phone

With the screen facing the lens, insert the bottom of the VR phone into the Device Holder. **Push the phone firmly into the plug** so you can't see the plug.



### 9. Lock phone in

Push down on the top of the phone (near the camera) and listen for the click to make sure the phone is locked in.



To remove phone, pull back on this latch and the phone will pop out.

## STEP 5: CONNECT THE HEARTRATE SENSOR

Wear the Scosche Rhythm+ sensor when you use the VR headset. The VR content will change based on your heart rate.

**Please note:** if your VR phone is not on, your sensor will not connect.

### TURN ON

Press and hold the button in the middle of the sensor for 2 seconds. The light next to the button will flash red and blue as it connects to the VR phone. Once it's connected, it will only flash blue.

If not properly connected, it'll flash red.

*Sensor Button*



### WEAR THE SENSOR

Included in your kit are 2 bands to fit your forearm.

Wear the armband on the upper part of your forearm. The sensor should be on the inside or top of your forearm. It should be tight enough not to move on your arm, but not so tight it restricts circulation.



# STEP 6: USE THE IPHONE APP

The EaseVR app on your iPhone will help guide you through the 10-day program.

## DOWNLOAD THE APP

Go to the App Store on your phone and search for “EaseVR” to download the app.



## FIRST-TIME USE

Open the EaseVR app on your iPhone and tap the Login link in the upper right corner.



**EMAIL:**

---

**PASSWORD:**


---



*Login*

## SURVEYS

Complete the questions until you reach the Home screen.

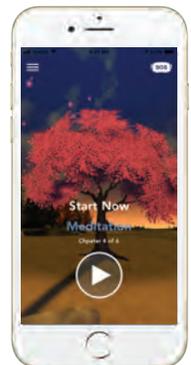

## STEP 7: START YOUR DAILY PRACTICE

Do you have the heartrate sensor on?

Are your headphones on?

Put the headphones around your neck so they are easy to access when you have the VR headset on. Make sure the microphone will be in front of your mouth.



Is the VR device ready?

Make sure you opened the PlayVR app on the VR phone and inserted the phone into the headset.



**Try to use the device in a quiet area with little background noise.**

1

**Press Play**

On your iPhone, press play for your daily chapter.



2

**Put on headset and headphones**

During the 10-second countdown, put on the VR headset. Then put on the headphones.

3

**Enjoy!**



Sit back, wait 10 seconds for it to start, and be immersed in your daily VR chapter.

## STEP 8: GET RELIEF

Use **Relief** when you have a pain flare and need immediate help for your pain symptoms.

Tap the Relief symbol in the upper right corner of the Home screen.

Pick from one of the individual sessions, tap Play and then put your VR headset on.

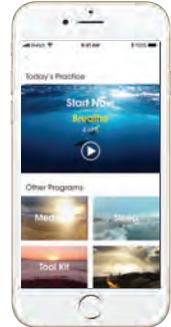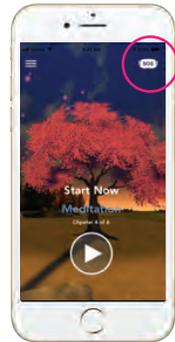


## STEP 9: TROUBLESHOOTING

### IPHONE APP AND VR DEVICE ARE NOT PAIRED?

If you press Play on your **iPhone app** and get an error message, please remove the **VR phone** from the headset.

Select the VR chapter on the VR phone.

Insert the VR phone into the headset and enjoy the VR session.

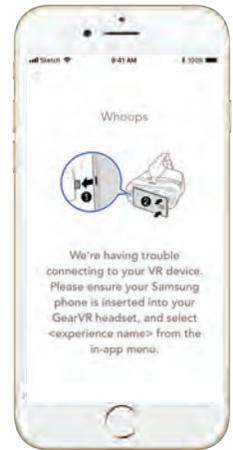

### HELPFUL TIPS


- Make sure your phone is connected to WIFI.
- Stay seated when using the VR headset. The VR headset is not intended to be used when walking around.
- If it's blurry, move the headset up and down your face to find the right focal point in the lense for a clearer view.
- Wear glasses inside the headset to improve your view.
- Place the microphone next to your mouth to ensure it picks up your breathing.
- Keep the VR device charged as the battery is continuously used when the VR phone is plugged into the headset.

# STEP 10: TURN OFF AND CHARGING

## HEART RATE SENSOR

Hold the button on the sensor down for 3 seconds. The light on the sensor will turn off.

To charge the sensor, connect the charging base to the back of the sensor. Plug the USB cord into the phone adapter. Charge for 2 hours or until a solid blue light appears on the sensor to indicate it is fully charged. The battery typically lasts for 8 hours.

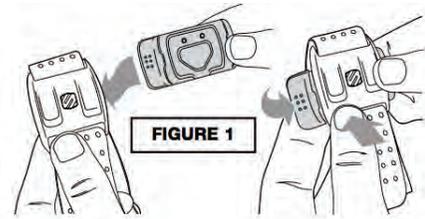

## CHARGE THE VR HEADSET

Your VR device will automatically go to sleep after a few minutes of inactivity. Plug the charger into the bottom of the headset to ensure your device has enough battery when you need to use it.

**If you do not charge the device, please remove the phone from the headset so the battery does not drain.**



*Plug the cord into the port on the underside of the headset*



# CONTACT US



[service@appliedvr.io](mailto:service@appliedvr.io)



844-857-0010



<https://appliedvr.zendesk.com/hc/en-us>